

PCT

WORLD INTELLECTUAL PROPERTY ORGANIZATION  
International Bureau



INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

<p>(51) International Patent Classification <sup>6</sup> : <b>G01N 33/53, 33/534</b></p>		<p>A1</p>	<p>(11) International Publication Number: <b>WO 97/10502</b> (43) International Publication Date: <b>20 March 1997 (20.03.97)</b></p>
<p>(21) International Application Number: <b>PCT/US96/14563</b> (22) International Filing Date: <b>11 September 1996 (11.09.96)</b></p>		<p>(81) Designated States: CA, JP, US, European patent (AT, BE, CH, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE).</p>	
<p>(30) Priority Data: 60/003,824 15 September 1995 (15.09.95) US 9603486.3 20 February 1996 (20.02.96) GB</p>		<p><b>Published</b> <i>With international search report. Before the expiration of the time limit for amending the claims and to be republished in the event of the receipt of amendments.</i></p>	
<p>(71) Applicant (<i>for all designated States except US</i>): MERCK &amp; CO., INC. [US/US]; 126 East Lincoln Avenue, Rahway, NJ 07065 (US).</p>			
<p>(72) Inventor; and (75) Inventor/Applicant (<i>for US only</i>): SALOWE, Scott, P. [US/US]; 126 East Lincoln Avenue, Rahway, NJ 07065 (US).</p>			
<p>(74) Common Representative: MERCK &amp; CO., INC.; 126 East Lincoln Avenue, Rahway, NJ 07065 (US).</p>			
<p>(54) Title: <b>A HIGH THROUGHPUT ASSAY USING FUSION PROTEINS</b></p>			
<p>(57) Abstract</p> <p>This application describes a high throughput assay for screening for compounds which are capable of binding to a fusion protein which consists of a target protein and an FK506-binding protein.</p>			

**FOR THE PURPOSES OF INFORMATION ONLY**

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AM	Armenia	GB	United Kingdom	MW	Malawi
AT	Austria	GE	Georgia	MX	Mexico
AU	Australia	GN	Guinea	NE	Niger
BB	Barbados	GR	Greece	NL	Netherlands
BE	Belgium	HU	Hungary	NO	Norway
BF	Burkina Faso	IE	Ireland	NZ	New Zealand
BG	Bulgaria	IT	Italy	PL	Poland
BJ	Benin	JP	Japan	PT	Portugal
BR	Brazil	KE	Kenya	RO	Romania
BY	Belarus	KG	Kyrgyzstan	RU	Russian Federation
CA	Canada	KP	Democratic People's Republic of Korea	SD	Sudan
CF	Central African Republic	KR	Republic of Korea	SE	Sweden
CG	Congo	KZ	Kazakhstan	SG	Singapore
CH	Switzerland	LI	Liechtenstein	SI	Slovenia
CI	Côte d'Ivoire	LK	Sri Lanka	SK	Slovakia
CM	Cameroon	LR	Liberia	SN	Senegal
CN	China	LT	Lithuania	SZ	Swaziland
CS	Czechoslovakia	LU	Luxembourg	TD	Chad
CZ	Czech Republic	LV	Latvia	TG	Togo
DE	Germany	MC	Monaco	TJ	Tajikistan
DK	Denmark	MD	Republic of Moldova	TT	Trinidad and Tobago
EE	Estonia	MG	Madagascar	UA	Ukraine
ES	Spain	ML	Mali	UG	Uganda
FI	Finland	MN	Mongolia	US	United States of America
FR	France	MR	Mauritania	UZ	Uzbekistan
GA	Gabon			VN	Viet Nam

TITLE OF THE INVENTION

## A HIGH THROUGHPUT ASSAY USING FUSION PROTEINS

BACKGROUND OF THE INVENTION

5        Src homology 2 (SH2) domains are a family of homologous protein domains that share the common property of recognizing phosphorylated tyrosine residues in specific peptide contexts. They have routinely been expressed in *E. coli* as fusion proteins with glutathione-S-transferase (GST). This usually provides high level expression and

10      straightforward affinity purification on glutathione-Sepharose. Ligand binding is then assayed by incubating the GST/S2 with a radiolabeled phosphopeptide, precipitating the complex with glutathione-Sepharose, washing the beads, and then counting the beads to determine bound radioactivity [Isakov et al., *J. Exp. Med.*, 181, 375-380 (1995); Piccione et al., *Biochemistry*, 32, 3197-3202 (1993); Huyer et al., *Biochemistry*, 34, 1040-1049 (1995)]. There are several disadvantages to this procedure, particularly when applied to high-throughput screening for agonists, antagonists, or inhibitors as new leads for drug development. First, the radiolabeling of the peptide is carried out either enzymatically

20      with a kinase and [<sup>32</sup>P]ATP or chemically with [<sup>125</sup>I]Bolton-Hunter reagent. In both cases, the isotopes are short-lived and thus require frequent preparation of material. In the case of enzymatic phosphorylation, the appropriate kinase must also be available in sufficient quantity to generate enough material for screening purposes.

25      Second, the protocol requires separation of bound complex from free phosphopeptide by washing of the glutathione-Sepharose beads. This is a nonequilibrium procedure that risks dissociation of the bound ligand, particularly when off-rates are fast. Thus, there is the possibility of misleading results. Finally, due to the number of manipulations and

30      centrifugations involved, the protocol is very tedious to conduct manually and is not readily adaptable to robotic automation to increase throughput.

Two additional methods for measuring the interaction of proteins and ligands that have been applied to SH2 domains are biospecific interaction analysis using surface plasmon resonance and

- 2 -

isothermal titration calorimetry (Felder et al., *Mol. Cell. Biol.*, 13, 1449-1455 (1993); Panayotou et al., *Mol. Cell. Biol.*, 13, 3567-3576 (1993); Payne et al., *Proc. Natl. Acad. Sci. U.S.A.*, 90, 4902-4906 (1993); Morelock et al., *J. Med. Chem.* 38, 1309-18 (1995); Ladbury et al., *Proc. Natl. Acad. Sci. U.S.A.*, 92, 3199-3203 (1995); Lemmon et al., *Biochemistry*, 33, 5070-5076 (1994)). These techniques do not require a particular fusion partner for the SH2 domain, but do require sophisticated instrumentation that is not amenable to high throughput screening.

10 **SUMMARY OF THE INVENTION**

The instant invention covers a method of screening for compounds capable of binding to a fusion protein which comprises combining a test compound, a tagged ligand, a fusion protein (target protein, peptide linker and FK506-binding protein), a radiolabeled ligand, and coated scintillation proximity assay (SPA) beads, and then measuring the scintillation counts attributable to the binding of the tagged ligand to the fusion protein in the presence of the test compound relative to a control assay in the absence of the test compound, so as to determine the effect the test compound has on the binding of the tagged ligand. This invention provides an immediate means of making use of SPA technology for the functional assay of ligand binding to a single or multiple signal transduction domain(s), for example a phosphopeptide binding to an SH2 domain. The present invention does not require specialized radiochemical synthesis and is readily adaptable to robotic automation for high capacity screening for agonists, antagonists, and/or inhibitors.

- 3 -

### BRIEF DESCRIPTION OF THE FIGURES

#### Figure 1.

5      A.) Binding of the streptavidin SPA bead, biotinylated ligand and the fusion protein (SH2:FKBP), which emits a detectable signal; and

5      B.) Binding of the test compound and the fusion protein (SH2:FKBP), which results in no signal detection .

### DETAILED DESCRIPTION OF THE INVENTION

10     The present invention relates to a method of screening for compounds which preferentially bind to a target protein.

An embodiment of this invention is a method of screening for compounds capable of binding to a fusion protein which comprises the steps of:

15     a) mixing a test compound, a tagged ligand, the fusion protein, a radiolabeled ligand and coated scintillation proximity assay (SPA) beads;

20     b) incubating the mixture for between about 1 hour and about 24 hours;

20     c) measuring the SPA bead-bound counts attributable to the binding of the tagged ligand to the fusion protein in the presence of the test compound using scintillation counting; and

25     d) determining the binding of the tagged ligand to the fusion protein in the presence of the test compound relative to a control assay run in the absence of the test compound.

The term "fusion protein" refers to a "target protein" fused to an "FK506-binding protein" (FKBP), the two proteins being separated by a "peptide linker".

30     A "peptide linker" may consist of a sequence containing from about 1 to about 20 amino acids, which may or may not include the sequence for a protease cleavage site. An example of a peptide linker which is a protease cleavage site is represented by the amino acid sequence GLPRGS.

- 4 -

The term "target protein" refers to any protein that has a defined ligand. Included within this definition of target protein are single and multiple signal transduction domains, such as, but not limited to, Src homology 1 (SH1), Src homology 2 (SH2), Src homology 3 (SH3), and 5 pleckstrin homology (PH) domains [Hanks & Hunter, *FASEB J.*, 9, 576-596 (1995); Bolen, *Curr. Opin. Immunol.*, 7, 306-311 (1995); Kuriyan & Cowburn, *Curr. Opin. Struct. Biol.*, 3, 828-837 (1993); Cohen et al., *Cell*, 80, 237-248 (1995)]. The term "SH1 domain" refers to a family of homologous protein domains that bind ATP and catalyze tyrosine 10 phosphorylation of peptide and protein substrates. The term "SH2 domain" refers to a family of homologous protein domains that share the common property of recognizing phosphorylated tyrosine residues in specific peptide contexts. The term "SH3 domain" refers to a family of homologous protein domains that share the common property of 15 recognizing polyproline type II helices. The term "PH domain" refers to a family of homologous protein domains that mediate both protein-protein and protein-lipid interactions. Examples of SH2 domains which may be utilized in the method of the invention include, but are not limited to, the single and tandem SH2 domains present in the tyrosine kinases 20 ZAP, SYK and LCK. The DNA sequences were obtained from GenBank, National Center for Biotechnology Information, National Library of Medicine, 8600 Rockville Pike, Bethesda, MD 20894. The Accession Numbers for the sequences are: human ZAP (L05148); human SYK (L28824) and human LCK (X13529).

25 The term "tagged ligand" refers to a biotinylated or epitope tagged ligand for the target protein.

The term "radiolabeled ligand" refers to a [<sup>3</sup>H]- or [<sup>125</sup>I]-labeled ligand which binds to the FKBP. An example of a radiolabeled ligand useful in the instant invention is [<sup>3</sup>H]-dihydroFK506.

30 The term "coated scintillation proximity assay beads" (SPA beads) refers to streptavidin-coated scintillation proximity assay beads when the tagged ligand is biotinylated, and to anti-epitope antibody bound to anti-antibody-coated or protein A-coated scintillation proximity assay beads when the tagged ligand is epitope-tagged.

- 5 -

The term "control assay" refers to the assay when performed in the presence of the tagged ligand, the fusion protein, the radiolabeled ligand and the coated scintillation proximity assay beads, but in the absence of the test compound.

5 The term FK506-binding proteins may include, but are not limited to, the below listed FKBP and FKBP homologues, which include a citation to the references which disclose them. This list is not intended to limit the scope of the invention.

10 **Mammalian**

FKBP-12 Galat et al., *Eur. J. Biochem.*, 216:689-707 (1993).  
FKBP-12.6 Wiederrecht, G. and F. Etzkorn  
*Perspectives in Drug Discovery and Design*, 2:57-84 (1994).  
15 FKBP-13 Galat et al., *supra*; Wiederrecht and Etzkorn, *supra*.  
FKBP-25 Galat et al., *supra*; Wiederrecht and Etzkorn, *supra*.  
20 FKBP-38 Wiederrecht and Etzkorn, *supra*.  
FKBP-51 Baughman et al., *Mol. Cell. Biol.*, 8, 4395-4402(1995).  
FKBP-52 Galat et al., *supra*.

25 **Bacteria**

Legionella pneumophilia Galat et al., *supra*.  
Legionella micadei Galat et al., *supra*.  
Chlamydia trachomatis Galat et al., *supra*.  
E. coli fkpa Horne, S.M. and K.D. Young, *Arch. Microbiol.*, 163:357-365 (1995).  
30 E. coli slyD Roof et al., *J. Biol. Chem.* 269:2902-2910 (1994).  
E. coli orf149 Trandinh et al., *FASEB J.* 6:3410-3420 (1992).

- 6 -

	Neisseria meningitidis	Hacker, J. and G. Fischer, <i>Mol. Micro.</i> , 10:445-456 (1993).
	Streptomyces chrysomallus	Hacker and Fischer, <i>supra</i> .
5	<b>Fungal</b>	
	yeast FKBP-12	Cardenas et al., <i>Perspectives in Drug Discovery and Design</i> , 2:103-126 (1994).
10	yeast FKBP-13	Cardenas et al., <i>supra</i> .
	yeast NPR1(FPR3)	Cardenas et al., <i>supra</i> .
	Neurospora	Galat et al., <i>supra</i> .

15 A variety of host cells may be used in this invention, which include, but are not limited to, bacteria, yeast, bluegreen algae, plant cells, insect cells and animal cells.

20 Expression vectors are defined herein as DNA sequences that are required for the transcription of cloned copies of genes and the translation of their mRNAs in an appropriate host. Such vectors can be used to express genes in a variety of host cells, such as, bacteria, yeast, bluegreen algae, plant cells, insect cells and animal cells.

25 Specifically designed vectors allow the shuttling of DNA between hosts such as bacteria-yeast or bacteria-animal cells. An appropriately constructed expression vector may contain: an origin of replication for autonomous replication in host cells, selectable markers, a limited number of useful restriction enzyme sites, a potential for high copy number, and active promoters. A promoter is defined as a DNA sequence that directs RNA polymerase to bind to DNA and initiate RNA synthesis. A strong promoter is one which causes mRNAs to be initiated at high frequency. Expression vectors may include, but are not limited 30 to, cloning vectors, modified cloning vectors, specifically designed plasmids or viruses. Commercially available vectors suitable for FKBP fusion protein expression include, but are not limited to pBR322 (Promega), pGEX (Amersham), pT7 (USB), pET (Novagen), pIB1 (IBI), pProEX-1 (Gibco/BRL), pBluescript II (Stratagene), pTZ18R and

- 7 -

pTZ19R (USB), pSE420 (Invitrogen), pVL1392 (Invitrogen), pBlueBac (Invitrogen), pBAcPAK (Clontech), pHIL (Invitrogen), pYES2 (Invitrogen), pCDNA (Invitrogen), pREP (Invitrogen) or the like.

The expression vector may be introduced into host cells via

- 5 any one of a number of techniques including but not limited to transformation, transfection, infection, protoplast fusion, and electroporation.

*E. coli* containing an expression plasmid with the target gene fused to FKBP are grown and appropriately induced. The cells are then

- 10 pelleted and resuspended in a suitable buffer. Although FKBP-12 lacks sequences that specifically direct it to the periplasm, FKBP fusions are primarily located there and can be released by a standard freeze/thaw treatment of the cell pellet. Following centrifugation, the resulting supernatant contains >80% pure FKBP fusion, which if desired can be
- 15 purified further by conventional methods. Alternatively, the assay is not dependent on pure protein and the initial periplasmic preparation may be used directly. A thrombin site located between FKBP and the target protein can be used as a means to cleave FKBP from the fusion; such cleaved material may be a suitable negative control for subsequent
- 20 assays.

A fusion protein which contains a single or multiple SH2 domain(s) may be purified by preparing an affinity matrix consisting of biotinylated phosphopeptide coupled to avidin or streptavidin immobilized on a solid support. A freeze/thaw extract is prepared from

- 25 the cells which express the fusion protein and is loaded onto the affinity matrix. The desired fusion protein is then specifically eluted with phenyl phosphate.

To assay the formation of a complex between a target protein and its ligand, the tagged ligand is mixed with the FKBP fusion protein in a suitable buffer in the presence of the radiolabeled ligand in the well of a white microplate. After a suitable incubation period to allow complex formation to occur, coated SPA beads are added to capture the tagged ligand and any bound fusion protein. The plate is sealed, incubated for a sufficient period to allow the capture to go to

- 30

- 8 -

completion, then counted in a multiwell scintillation counter. Screening for agonists/antagonists/inhibitors is carried out by performing the initial incubation prior to the capture step with SPA beads in the presence of a test compound(s) to determine whether they have an effect upon the  
5 binding of the tagged ligand to the fusion protein. This principle is illustrated by Figure 1.

The present invention can be understood further by the following examples, which do not constitute a limitation of the invention.

10

### EXAMPLE 1

#### Process for Preparing the FKBP fusion cloning vector

General techniques for modifying and expressing genes in various host cells can be found in Ausubel, F.M., Brent, R., Kingston, R.E., Moore, D.D., Seidman, J.G., Smith, J.A. and Struhl, K. *Current Protocols in Molecular Biology* (John Wiley & Sons, New York, New York, 1989). Sequences for a 3'- altered FKBP fragment that contained a glycine codon (GGT) in place of the stop (TGA) codon followed by a sequence encoding a thrombin site (Leu-Val-Pro-Arg) and *Bam*HI  
15 restriction site (GAATTC) were amplified using the polymerase chain reaction (PCR). The PCR reaction contained the following primers: 5'- GATGCCATGGGAGTGCAGGTGGAAACCATCTCCCCA-3' and 5'- TACGAATTCTGGCGTGGATCCACGCCGAACCAGACCTTCCAGT TTTAG-3' and a plasmid containing human FKBP-12 as the template.  
20 The resulting 367 base pair amplification product was ligated into the vector pCRII (Invitrogen) and the ligation mixture transformed into competent *Escherichia coli* cells. Clones containing an insert were identified using PCR with flanking vector primers. Dideoxy DNA sequencing confirmed the nucleotide sequence of one positive isolate.  
25 The altered 338 base pair FKBP fragment was excised from the pCRII plasmid using *Nco*I and *Bam*HI and ligated into *Nco*I and *Bam*HI digested pET9d (Novagen) plasmid. Competent *E. coli* were transformed with the ligation mixture, and colonies containing the insert were identified using PCR with primers encoding for flanking vector sequences. The FKBP  
30 fusion cloning vector is called pET9dFKBPt.  
35

EXAMPLE 2Process for Preparing the FK-ZAP fusion expression vector

5 A DNA fragment encoding for the tandem SH2 domains of ZAP-70 was prepared by PCR to contain a *Bam*HI site at the 5'-end such that the reading frame was conserved with that of FKBP in the fusion vector. At the 3'-end, the fragment also incorporated a stop codon followed by a *Bam*HI site. The PCR reaction contained Molt-4 cDNA  
10 (Clontech) and the following primers:  
5'-ATTAGGATCCATGCCAGATCCTGCAGCTCACCTGCCCT-3' and  
5'-ATATGGATCCTTACCAGAGGCGTTGCT-3'. The fragment was cloned into a suitable vector, sequenced, digested with *Bam*HI, and the insert containing the SH2 domains ligated to *Bam*HI treated  
15 pET9dFKB Pt, and transformed into *E. coli*. Clones containing inserts in the correct orientation were identified by PCR or restriction enzyme analysis. Plasmid DNA was prepared and used to transform BL21(DE3) cells.

20 EXAMPLE 3Process for Preparing the FK-SYK fusion expression vector

The expression vector for the tandem SH2 domains of Syk fused to FKBP was prepared as in Example 2 except that the PCR  
25 reaction contained Raji cell cDNA (Clontech) and the following primers:  
5'-CAATAGGATCCATGCCAGCAGCGGCATGGCTGA-3' and  
5'-GACCTAGGATCCCTAATTAACATTTCCCTGTGTGCCGAT-3'.

EXAMPLE 4

30

Process for Preparing the FK-LCK fusion expression vector

The expression vector for the SH2 domain of Lck fused to FKBP was prepared as in Example 2 except that the PCR reaction contained Molt-4 cDNA (Clontech) and the following primers:

- 10 -

5'-ATATGGATCCATGGCGAACAGCCTGGAGCCCGAACCT-3'  
and  
5'-ATTAGGATCCTTAGGTCTGGCAGGGCGGCTAACCGTGT  
GCA-3'.

5

#### EXAMPLE 5

##### FK-ZAP

10 Step A: Process for Expression of FK-ZAP  
*E. coli* BL21(DE3) cells containing the pET9dFKB Pt/  
ZapSH2 plasmid were grown in Luria-Bertani (LB) media containing 50  
microgram/ml kanamycin at 37 degrees C until the optical density  
measured at 600 nm was 0.5-1.0. Expression of the FK-ZAP fusion  
15 protein was induced with 0.1 mM isopropyl beta-thiogalactopyranoside  
and the cells were grown for another 3-5 hr at 30 degrees C. They were  
pelleted at 4400 x g for 10 min at 4 degrees C and resuspended in 2% of  
the original culture volume with 100 mM tris pH 8.0 containing 1  
microgram/ml each aprotinin, pepstatin, leupeptin, and bestatin. The  
20 resuspended pellet was frozen at -20 degrees C until further purification.

Step B: Process for Purification of FK-ZAP  
The affinity matrix for purification of FK-ZAP was prepared  
by combining agarose-immobilized avidin with excess biotinylated  
25 phosphopeptide derived from the  $\zeta 1$  ITAM sequence of the human T-cell  
receptor, biotinyl-GSNQLpYNELNLGRREEpYDVLDK, and washing  
out unbound peptide. Frozen cells containing FK-ZAP were thawed in  
warm water, refrozen on dry ice for about 25 min., then thawed again.  
After the addition of 0.1% octyl glucoside, 1 mM dithiothreitol (DTT)  
30 and 500 mM NaCl, the extract was centrifuged at 35,000 x g for  
approximately 30 minutes. The supernatant was loaded onto the  
phosphopeptide affinity column, at about 4° and washed with phosphate  
buffered saline containing 1 mM DTT and 0.1% octyl glucoside.  
FK-ZAP was eluted with 200 mM phenyl phosphate in the same buffer at  
35 about 37°. The protein pool was concentrated and the phenyl phosphate

- 11 -

removed on a desalting column. The purified FK-ZAP was stored at about -30° in 10 mM HEPES/150 mM NaCl/1 mM DTT/0.1 mM EDTA/10% glycerol.

5

EXAMPLE 6

FK-SYK

10 *E. coli* BL21(DE3) cells containing the pET9dFKB Pt/  
SykSH2 plasmid were grown, induced, and harvested as described in  
Example 5. FK-SYK was purified using the same affinity matrix and  
methodology described in Example 5.

EXAMPLE 7

15 FK-LCK

20 *E. coli* BL21(DE3) cells containing the pET9dFKB Pt/  
LckSH2 plasmid were grown, induced, and harvested as described in  
Example 5. The affinity matrix for purification of FK-LCK was prepared  
by combining agarose-immobilized avidin with excess biotinyl-  
EPQpYEEIPIYL, and washing out unbound peptide. The remaining  
methodology for purification was the same as Example 5.

EXAMPLE 8

25

Method of Screening for Antagonists of FK-ZAP

30 Assays were conducted at ambient temperature in a buffer  
consisting of 25 mM HEPES, 10 mM DTT, 0.01% TWEEN-20, pH 7.0.  
10 µl of a DMSO solution of test compound(s) and 120 µl of biotinyl-  
phosphopeptide stock solution were dispensed into the wells of a 96-well  
Packard Optiplate. Next, 20 µl of a mixture of FK-ZAP protein and  
<sup>3</sup>H-dihydroFK506 were added to each test well. Finally, 50 µl of a 4  
mg/ml suspension of SPA beads were dispensed to each well. Final  
concentrations of the assay components were:

- 12 -

25 nM biotinyl-GSNQLpYNELNLGRREEpYDVLDK  
25 nM FK-ZAP fusion protein  
10 nM  $^3$ H-dihydroFK506 (DuPont NEN)  
1.0 mg/ml streptavidin-SPA beads (Amersham)  
5% DMSO

5 The plate was sealed and incubated between 1 and 8 hours. Bead-bound radioactivity was then measured in a Packard Topcount microplate scintillation counter.

10 EXAMPLE 9

Method of Screening for Antagonists of FK-SYK

The assays were conducted as set forth in Example 8, except that FK-SYK replaced FK-ZAP.

15 EXAMPLE 10

Method of Screening for Antagonists of FK-LCK

20 The assays were conducted as set forth in Example 8, except that FK-LCK replaced FK-ZAP and the tagged ligand was 25 nM biotinyl-EPQpYEEIPIYL.

- 13 -

SEQUENCE LISTING

(1) GENERAL INFORMATION:

- (i) APPLICANT: Salowe, Scott P.
- (ii) TITLE OF INVENTION: A High Throughput Assay Using Fusion Proteins
- (iii) NUMBER OF SEQUENCES: 6
- (iv) CORRESPONDENCE ADDRESS:
  - (A) ADDRESSEE: Valerie J. Camara
  - (B) STREET: 126 E. Lincoln Avenue, P.O. Box 2000
  - (C) CITY: Rahway
  - (D) STATE: NJ
  - (E) COUNTRY: U.S.A.
  - (F) ZIP: 07065
- (v) COMPUTER READABLE FORM:
  - (A) MEDIUM TYPE: Floppy disk
  - (B) COMPUTER: IBM PC compatible
  - (C) OPERATING SYSTEM: PC-DOS/MS-DOS
  - (D) SOFTWARE: PatentIn Release #1.0, Version #1.30
- (vi) CURRENT APPLICATION DATA:
  - (A) APPLICATION NUMBER:
  - (B) FILING DATE:
  - (C) CLASSIFICATION:
- (viii) ATTORNEY/AGENT INFORMATION:
  - (A) NAME: Camara, Valerie J.
  - (B) REGISTRATION NUMBER: 35,090
  - (C) REFERENCE/DOCKET NUMBER: 19494
- (ix) TELECOMMUNICATION INFORMATION:
  - (A) TELEPHONE: (908) 594-3902
  - (B) TELEFAX: (908) 594-4720

(2) INFORMATION FOR SEQ ID NO:1:

- (i) SEQUENCE CHARACTERISTICS:
  - (A) LENGTH: 1137 base pairs
  - (B) TYPE: nucleic acid
  - (C) STRANDEDNESS: single
  - (D) TOPOLOGY: linear
- (ii) MOLECULE TYPE: DNA (genomic)
- (xi) SEQUENCE DESCRIPTION: SEQ ID NO:1:

- 14 -

ATGGGAGTGC AGGTGGAAAC CATCTCCCCA GGAGATGGAC GCACCTTCCC CAAGCGCGC	60
CAGACCTGCG TGGTGCCTA CACCGGGATG CTTGAAGATG GAAAGAAATT TGATTCCTCC	120
CGGGACAGAA ACAAGCCCTT TAAGTTATG CTAGGCAAGC AGGAGGTGAT CCGAGGCTGG	180
GAAGAAGGGG TTGCCCAGAT GAGTGTGGGT CAGAGGCCA AACTGACTAT ATCTCCAGAT	240
TATGCCTATG GTGCCACTGG GCACCCAGGC ATCATCCAC CACATGCCAC TCTCGTCTTC	300
GATGTGGAGC TTCTAAAATC GGAAGGTCTG GTTCCGCGTG GATCCATGCC AGATCCTGCA	360
GCTCACCTGC CCTTCTTCTA CGGCAGCAGTC TCGCGTGCCTG AGGCCGAGGA GCACCTGAAG	420
CTGGCGGGCA TGGCGGACGG GCTCTTCTG CTGCGCCAGT GCCTGCGCTC GCTGGCGGGC	480
TATGTGCTGT CGCTCGTGCA CGATGTGCGC TTCCACCACT TTCCCATCGA GCGCCAGCTC	540
AACGGCACCT ACGCCATTGC CGGGGGCAAA GCGCACTGTG GACCGGCAGA GCTCTGCGAG	600
TTCTACTCGC GCGACCCCGA CGGGCTGCCCT TGCAACCTGC GCAAGCCGTG CAACCGGCCG	660
TCGGGCCTCG AGCCGCAGCC GGGGGTCTTC GACTGCCTGC GAGACGCCAT GGTGCGTGAC	720
TACGTGCGCC AGACGTGGAA GCTGGAGGGC GAGGCCCTGG AGCAGGCCAT CATCAGCCAG	780
GCCCCGCAGG TGGAGAAGCT CATTGCTACG ACGGCCACG AGCGGATGCC CTGGTACCAAC	840
AGCAGCCTGA CGCGTGAGGA GGCGAGCGT AAACTTTACT CTGGGGCGCA GACCGACGGC	900
AAGTTCTGC TGAGGCCCG GAGGAGCAG GGCACATACG CCCTGTCCCT CATCTATGGG	960
AAGACGGTGT ACCACTACCT CATCAGCCAA GACAAGGCCG GCAAGTACTG CATTCCGAG	1020
GGCACCAAGT TTGACACGCT CTGGCAGCTG GTGGAGTATC TGAAGCTGAA GGCGGACGGG	1080
CTCATCTACT GCCTGAAGGA GGCTGCCCT AACAGCAGTG CCAGCAACGC CTCTTAA	1137

**(2) INFORMATION FOR SEQ ID NO:2:**

- (i) SEQUENCE CHARACTERISTICS:**
  - (A) LENGTH: 1155 base pairs
  - (B) TYPE: nucleic acid
  - (C) STRANDEDNESS: single
  - (D) TOPOLOGY: linear

**(ii) MOLECULE TYPE: DNA (genomic)**

**(xi) SEQUENCE DESCRIPTION: SEQ ID NO:2:**

ATGGGAGTGC AGGTGGAAAC CATCTCCCCA GGAGATGGAC GCACCTTCCC CAAGCGCGC	60
CAGACCTGCG TGGTGCCTA CACCGGGATG CTTGAAGATG GAAAGAAATT TGATTCCTCC	120

- 15 -

CGGGACAGAA ACAAGCCCTT TAAGTTTATG CTAGGCAAGC AGGAGGTGAT CCGAGGCTGG	180
GAAGAAGGGG TTGCCCAGAT GAGTGTGGGT CAGAGAGCCA AACTGACTAT ATCTCCAGAT	240
TATGCCCTATG GTGCCACTGG GCACCCAGGC ATCATCCAC CACATGCCAC TCTCGTCTTC	300
GATGTGGAGC TTCTAAAACG GGAAGGTCTG GTTCCGCGTG GATCCATGGC CAGCAGCGGC	360
ATGGCTGACA GCGCCAACCA CCTGCCCTTC TTTTTCGGCA ACATCACCCG GGAGGAGGCA	420
GAAGATTACC TGGTCCAGGG GGGCATGAGT GATGGGCTTT ATTTGCTGCG CCAGAGCCGC	480
AACTACCTGG GTGGCTTCGG CCTGTCCGTG GCCCACGGGA GGAAGGCACA CCACTACACCC	540
ATCGAGCGGG AGCTGAATGG CACCTACGCC ATCGCCGGTG GCAGGACCA TGCCAGCCCC	600
GCCGACCTCT GCCACTACCA CTCCCAGGAG TCTGATGGCC TGGTCTGCCT CCTCAAGAAG	660
CCCTTCAACC GGCCCCAAGG GGTGCAGCCC AAGACTGGC CCTTTGAGGA TTTGAAGGAA	720
AACCTCATCA GGGAAATATGT GAAGCAGACA TGGAACCTGC AGGGTCAGGC TCTGGAGCAG	780
GCCATCATCA GTCAGAAGCC TCAGCTGGAG AAGCTGATCG CTACCACAGC CCATGAAAAA	840
ATGCCTTGGT TCCATGGAAA AATCTCTCGG GAAGAATCTG AGCAAATTGT CCTGATAGGA	900
TCAAAGACAA ATGGAAAGTT CCTGATCCGA GCCAGAGACA ACAACGGCTC CTACGCCCTG	960
TGCCTGCTGC ACGAAGGGAA GGTGCTGCAC TATCGCATCG ACAAAAGACAA GACAGGGAAAG	1020
CTCTCCATCC CCGAGGGAAA GAAGTTCGAC ACCGCTCTGGC AGCTAGTCGA GCATTATTCT	1080
TATAAAGCAG ATGGTTTGGT AAGAGTTCTT ACTGTCCCAT GTCAAAAAT CGGCACACAG	1140
GGAAATGTTA ATTAG	1155

## (2) INFORMATION FOR SEQ ID NO:3:

(i) SEQUENCE CHARACTERISTICS:

- (A) LENGTH: 675 base pairs
- (B) TYPE: nucleic acid
- (C) STRANDEDNESS: single
- (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: DNA (genomic)

## (xi) SEQUENCE DESCRIPTION: SEQ ID NO:3:

ATGGGAGTGC AGGTGGAAAC CATCTCCCCA GGAGATGGAC GCACCTTCCC CAAGCGCGGC	60
CAGACCTGCG TGGTGCACTA CACCGGGATG CTTGAAGATG GAAAGAAATT TGATTCTCC	120

- 16 -

CGGGACAGAA ACAAGCCCTT TAAGTTTATG CTAGGCAAGC AGGAGGTGAT CCGAGGCTGG	180
GAAGAAGGGG TTGCCCCAGAT GAGTGTGGGT CAGAGAGCCA AACTGACTAT ATCTCCAGAT	240
TATGCCTATG GTGCCACTGG GCACCCAGGC ATCATCCCAC CACATGCCAC TCTCGTCTTC	300
GATGTGGAGC TTCTAAAATC GGAAGGTCTG GTTCCGCGTG GATCCATGGC GAACAGCCTG	360
GAGCCCGAAC CCTGGTTCTT CAAGAACCTG AGCCGCAAGG ACGCGGAGCG GCAGCTCCTG	420
GCGCCCGGGA ACACTCACGG CTCCCTTCCTC ATCCGGGAGA GCGAGAGCAC CGCGGGATCG	480
TTTTCACTGT CGGTCCGGGA CTTCGACCAG AACCAAGGAG AGGTGGTGAA ACATTACAAG	540
ATCCGTAATC TGGACAACGG TGGCTTCTAC ATCTCCCTC GAATCACTTT TCCCGGCCTG	600
CATGAACTGG TCCGCCATTA CACCAATGCT TCAGATGGGC TGTGCACACG GTTGAGCCGC	660
CCCTGCCAGA CCTAA	675

(2) INFORMATION FOR SEQ ID NO:4:

**(i) SEQUENCE CHARACTERISTICS:**

- (A) LENGTH: 378 amino acids
- (B) TYPE: amino acid
- (C) STRANDEDNESS: single
- (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: protein

Met Gly Val Gln Val Glu Thr Ile Ser Pro

Pre-Law Law School Test Study Guide 15

20 25 30

Phe Met Leu Gly Lys Gln Glu Val Ile Arg Gly Trp Glu Glu Gly Val  
50 55 60

Ala Gln Met Ser Val Gly Gln Arg Ala Lys Leu Thr Ile Ser Pro Asp  
65 70 75 80

Tyr Ala Tyr Gly Ala Thr Gly His Pro Gly Ile Ile Pro Pro His Ala  
85 90 95

Thr Leu Val Phe Asp Val Glu Leu Leu Lys Leu Glu Gly Leu Val Pro  
100 105 110

- 17 -

Arg Gly Ser Met Pro Asp Pro Ala Ala His Leu Pro Phe Phe Tyr Gly  
115 120 125

Ser Ile Ser Arg Ala Glu Ala Glu Glu His Leu Lys Leu Ala Gly Met  
130 135 140

Ala Asp Gly Leu Phe Leu Leu Arg Gln Cys Leu Arg Ser Leu Gly Gly  
145 150 155 160

Tyr Val Leu Ser Leu Val His Asp Val Arg Phe His His Phe Pro Ile  
165 170 175

Glu Arg Gln Leu Asn Gly Thr Tyr Ala Ile Ala Gly Gly Lys Ala His  
180 185 190

Cys Gly Pro Ala Glu Leu Cys Glu Phe Tyr Ser Arg Asp Pro Asp Gly  
195 200 205

Leu Pro Cys Asn Leu Arg Lys Pro Cys Asn Arg Pro Ser Gly Leu Glu  
210 215 220

Pro Gln Pro Gly Val Phe Asp Cys Leu Arg Asp Ala Met Val Arg Asp  
225 230 235 240

Tyr Val Arg Gln Thr Trp Lys Leu Glu Gly Glu Ala Leu Glu Gln Ala  
245 250 255

Ile Ile Ser Gln Ala Pro Gln Val Glu Lys Leu Ile Ala Thr Thr Ala  
260 265 270

His Glu Arg Met Pro Trp Tyr His Ser Ser Leu Thr Arg Glu Glu Ala  
275 280 285

Glu Arg Lys Leu Tyr Ser Gly Ala Gln Thr Asp Gly Lys Phe Leu Leu  
290 295 300

Arg Pro Arg Lys Glu Gln Gly Thr Tyr Ala Leu Ser Leu Ile Tyr Gly  
305 310 315 320

Lys Thr Val Tyr His Tyr Leu Ile Ser Gln Asp Lys Ala Gly Lys Tyr  
325 330 335

Cys Ile Pro Glu Gly Thr Lys Phe Asp Thr Leu Trp Gln Leu Val Glu  
340 345 350

Tyr Leu Lys Leu Lys Ala Asp Gly Leu Ile Tyr Cys Leu Lys Glu Ala  
355 360 365

Cys Pro Asn Ser Ser Ala Ser Asn Ala Ser  
370 375

(2) INFORMATION FOR SEQ ID NO:5:

(i) SEQUENCE CHARACTERISTICS:

(A) LENGTH: 384 amino acids  
(B) TYPE: amino acid

- 18 -

(C) STRANDEDNESS: single  
 (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: protein

(xi) SEQUENCE DESCRIPTION: SEQ ID NO:5:

Met	Gly	Val	Gln	Val	Glu	Thr	Ile	Ser	Pro	Gly	Asp	Gly	Arg	Thr	Phe
1				5					10				15		
Pro	Lys	Arg	Gly	Gln	Thr	Cys	Val	Val	His	Tyr	Thr	Gly	Met	Leu	Glu
					20			25				30			
Asp	Gly	Lys	Lys	Phe	Asp	Ser	Ser	Arg	Asp	Arg	Asn	Lys	Pro	Phe	Lys
				35			40				45				
Phe	Met	Leu	Gly	Lys	Gln	Glu	Val	Ile	Arg	Gly	Trp	Glu	Glu	Gly	Val
		50				55				60					
Ala	Gln	Met	Ser	Val	Gly	Gln	Arg	Ala	Lys	Leu	Thr	Ile	Ser	Pro	Asp
		65			70			75			80				
Tyr	Ala	Tyr	Gly	Ala	Thr	Gly	His	Pro	Gly	Ile	Ile	Pro	Pro	His	Ala
				85			90				95				
Thr	Leu	Val	Phe	Asp	Val	Glu	Leu	Leu	Lys	Leu	Glu	Gly	Leu	Val	Pro
				100			105				110				
Arg	Gly	Ser	Met	Ala	Ser	Ser	Gly	Met	Ala	Asp	Ser	Ala	Asn	His	Leu
		115				120			125						
Pro	Phe	Phe	Phe	Gly	Asn	Ile	Thr	Arg	Glu	Glu	Ala	Glu	Asp	Tyr	Leu
	130					135			140						
Val	Gln	Gly	Gly	Met	Ser	Asp	Gly	Leu	Tyr	Leu	Leu	Arg	Gln	Ser	Arg
	145				150				155			160			
Asn	Tyr	Leu	Gly	Gly	Phe	Ala	Leu	Ser	Val	Ala	His	Gly	Arg	Lys	Ala
				165				170			175				
His	His	Tyr	Thr	Ile	Glu	Arg	Glu	Leu	Asn	Gly	Thr	Tyr	Ala	Ile	Ala
				180			185			190					
Gly	Gly	Arg	Thr	His	Ala	Ser	Pro	Ala	Asp	Leu	Cys	His	Tyr	His	Ser
		195				200				205					
Gln	Glu	Ser	Asp	Gly	Leu	Val	Cys	Leu	Leu	Lys	Lys	Pro	Phe	Asn	Arg
		210			215				220						
Pro	Gln	Gly	Val	Gln	Pro	Lys	Thr	Gly	Pro	Phe	Glu	Asp	Leu	Lys	Glu
	225				230			235			240				

- 19 -

Asn Leu Ile Arg Glu Tyr Val Lys Gln Thr Trp Asn Leu Gln Gly Gln  
245 250 255

Ala Leu Glu Gln Ala Ile Ile Ser Gln Lys Pro Gln Leu Glu Lys Leu  
260 265 270

Ile Ala Thr Thr Ala His Glu Lys Met Pro Trp Phe His Gly Lys Ile  
275 280 285

Ser Arg Glu Glu Ser Glu Gln Ile Val Leu Ile Gly Ser Lys Thr Asn  
290 295 300

Gly Lys Phe Leu Ile Arg Ala Arg Asp Asn Asn Gly Ser Tyr Ala Leu  
305 310 315 320

Cys Leu Leu His Glu Gly Lys Val Leu His Tyr Arg Ile Asp Lys Asp  
325 330 335

Lys Thr Gly Lys Leu Ser Ile Pro Glu Gly Lys Lys Phe Asp Thr Leu  
340 345 350

Trp Gln Leu Val Glu His Tyr Ser Tyr Lys Ala Asp Gly Leu Leu Arg  
355 360 365

Val Leu Thr Val Pro Cys Gln Lys Ile Gly Thr Gln Gly Asn Val Asn  
370 375 380

(2) INFORMATION FOR SEQ ID NO:6:

(i) SEQUENCE CHARACTERISTICS:

- (A) LENGTH: 224 amino acids
- (B) TYPE: amino acid
- (C) STRANDEDNESS: single
- (D) TOPOLOGY: linear

(ii) MOLECULE TYPE: protein

(xi) SEQUENCE DESCRIPTION: SEQ ID NO:6:

Met Gly Val Gln Val Glu Thr Ile Ser Pro Gly Asp Gly Arg Thr Phe  
1 5 10 15

Pro Lys Arg Gly Gln Thr Cys Val Val His Tyr Thr Gly Met Leu Glu  
20 25 30

Asp Gly Lys Lys Phe Asp Ser Ser Arg Asp Arg Asn Lys Pro Phe Lys  
35 40 45

Phe Met Leu Gly Lys Gln Glu Val Ile Arg Gly Trp Glu Glu Gly Val  
50 55 60

- 20 -

Ala Gln Met Ser Val Gly Gln Arg Ala Lys Leu Thr Ile Ser Pro Asp  
65 70 75 80

Tyr Ala Tyr Gly Ala Thr Gly His Pro Gly Ile Ile Pro Pro His Ala  
85 90 95

Thr Leu Val Phe Asp Val Glu Leu Leu Lys Leu Glu Gly Leu Val Pro  
100 105 110

Arg Gly Ser Met Ala Asn Ser Leu Glu Pro Glu Pro Trp Phe Phe Lys  
115 120 125

Asn Leu Ser Arg Lys Asp Ala Glu Arg Gln Leu Leu Ala Pro Gly Asn  
130 135 140

Thr His Gly Ser Phe Leu Ile Arg Glu Ser Glu Ser Thr Ala Gly Ser  
145 150 155 160

Phe Ser Leu Ser Val Arg Asp Phe Asp Gln Asn Gln Gly Glu Val Val  
165 170 175

Lys His Tyr Lys Ile Arg Asn Leu Asp Asn Gly Gly Phe Tyr Ile Ser  
180 185 190

Pro Arg Ile Thr Phe Pro Gly Leu His Glu Leu Val Arg His Tyr Thr  
195 200 205

Asn Ala Ser Asp Gly Leu Cys Thr Arg Leu Ser Arg Pro Cys Gln Thr  
210 215 220

WHAT IS CLAIMED IS:

1. A method of screening for compounds capable of binding to a fusion protein which comprises the steps of:
  - 5 a) mixing a test compound, a tagged ligand, the fusion protein, a radiolabeled ligand and coated scintillation proximity assay (SPA) beads;
  - b) incubating the mixture from between about 1 hour to about 24 hours;
  - 10 c) measuring the SPA bead-bound counts attributable to the binding of the tagged ligand to the fusion protein in the presence of the test compound using scintillation counting; and
  - d) determining the binding of the tagged ligand to the fusion protein in the presence of the test compound relative to a control assay run in the absence of the test compound.
- 20 2. The method of screening for compounds capable of binding to a fusion protein, as recited in Claim 1, wherein the tagged ligand is a biotinylated ligand or epitope-tagged ligand.
- 25 3. The method of screening for compounds capable of binding to a fusion protein, as recited in Claim 2, wherein scintillation proximity assay beads are streptavidin-coated or anti-antibody or protein A-coated.
- 30 4. The method of screening for compounds capable of binding to a fusion protein, as recited in Claim 3, wherein the radiolabeled ligand consists of [<sup>3</sup>H]- or [<sup>125</sup>I]-labeled FK506 analog.
5. The method of screening for compounds capable of binding to a fusion protein, as recited in Claim 4, wherein the fusion protein comprises an FK506-binding protein linked through a peptide linker to a target protein.

- 22 -

6. The method of screening for compounds capable of binding to a fusion protein, as recited in Claim 5, wherein the target protein comprises a single or multiple signal transduction domain.

5

7. The method of screening for compounds capable of binding to a fusion protein, as recited in Claim 6, wherein the single or multiple signal transduction domain is selected from the group consisting of: SH1, SH2, SH3 and PH domains.

10

8. The method of screening for compounds capable of binding to a fusion protein, as recited in Claim 7, wherein the target protein is a single or multiple SH2 domain.

15

9. The method of screening for compounds capable of binding to a fusion protein, as recited in Claim 8, wherein the radiolabeled ligand is [<sup>3</sup>H]-dihydroFK506.

20

10. The method of screening for compounds capable of binding to a fusion protein, as recited in Claim 9, wherein the FK506-binding protein is a 12kDA human FK506-binding protein.

25

11. The method of screening for compounds capable of binding to a fusion protein, as recited in Claim 10, wherein the target protein is a single or multiple SH2 domain selected from the group consisting of: ZAP:SH2, SYK:SH2 and LCK:SH2.

30

12. The method of screening for compounds capable of binding to a fusion protein, as recited in Claim 11, wherein the target protein is the SH2 domain, ZAP:SH2.

13. The method of screening for compounds capable of binding to a fusion protein, as recited in Claim 11, wherein the target protein is the SH2 domain, SYK:SH2.

- 23 -

14. The method of screening for compounds capable of binding to a fusion protein, as recited in Claim 11, wherein the target protein is the SH2 domain, LCK:SH2.

5

1/1

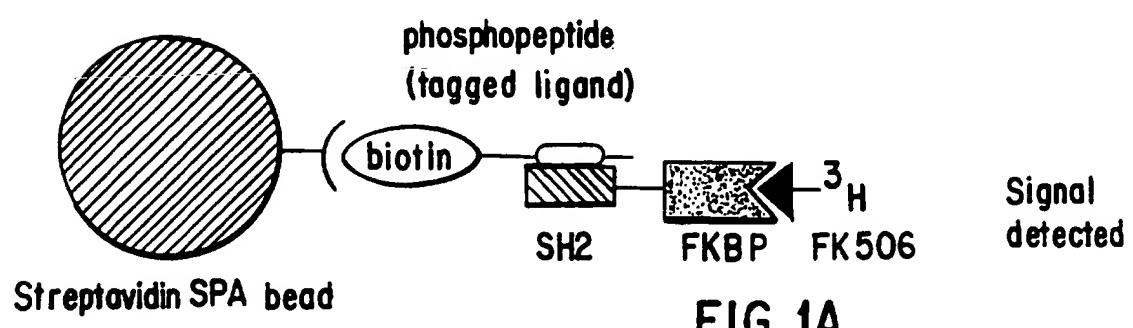


FIG. 1A

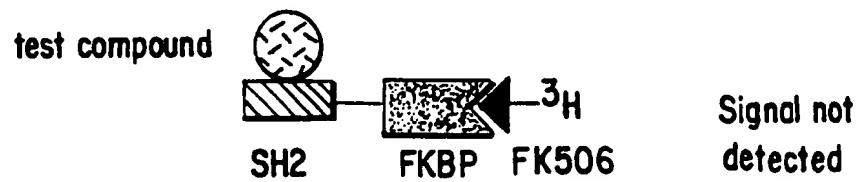


FIG. 1B

# INTERNATIONAL SEARCH REPORT

Int'l. application No.  
PCT/US96/14563

## A. CLASSIFICATION OF SUBJECT MATTER

IPC(6) : G01N 33/53, 33/534

US CL : 435/7.5; 436/529

According to International Patent Classification (IPC) or to both national classification and IPC

## B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)

U.S. : 435/7.5; 436/529

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practicable, search terms used)

APS, MEDLINE, EMBASE, BIOSIS, CAPLUS

search terms: scintillation proximity assay#, spa, fusion protein#, sh2 domain#, transduction domain#, target protein#, biotin?

## C. DOCUMENTS CONSIDERED TO BE RELEVANT

Category*	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X	Sonatore et al. The Utility of FK506-Binding Protein as a Fusion Partner in Scintillation Proximity Assays: Application to SH2 Domains. Anal. Biochem. 1996, Vol. 240, pages 289-297.	1-14

Further documents are listed in the continuation of Box C.  See patent family annex.

Special categories of cited documents:	
*A*	document defining the general state of the art which is not considered to be of particular relevance
*E*	earlier document published on or after the international filing date
*L*	document which may throw doubt on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)
*O*	document referring to an oral disclosure, use, exhibition or other means
*P*	document published prior to the international filing date but later than the priority date claimed
*T*	later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention
*X*	document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone
*Y*	document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art
*A*	document member of the same patent family

Date of the actual completion of the international search

08 DECEMBER 1996

Date of mailing of the international search report

10 JAN 1997

Name and mailing address of the ISA/US  
Commissioner of Patents and Trademarks  
Box PCT  
Washington, D.C. 20231

Faxsimile No. (703) 305-3230

Authorized officer

ROSEMARY ASHTON

Telephone No. (703) 308-0196

**INTERNATIONAL SEARCH REPORT**International application No.  
PCT/US96/14563**C (Continuation). DOCUMENTS CONSIDERED TO BE RELEVANT**

Category*	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X	Lerner et al. Scintillation Proximity Assay for Human DNA Topoisomerase I using Recombinant Biotiny Fusion Protein Produced in Baculovirus-Infected Insect Cells. <i>Anal. Biochem.</i> 1996, Vol. 240, pages 289-297.	
Y	Skinner et al. Direct Measurement of the Binding of RAS to Neurobromin using a Scintillation Proximity Assay. <i>Anal. Biochem.</i> 1994, Vol. 223, pages 259-265.	1,3